

IS-199 CHEMICAL MIXER



EXPLODED VIEW OF MIXER

FEATURES

Specific Gravity Mixing System

The Ultimate in Convenience
and Accuracy



Rugged Field Tested Components

Assure Reliability



Low Level Warning Systems

Audible and Visible



Easy-to-Service Design



UL/CSA Listed Power Supply

External Wall Mount
Class 2 Transformer




24 Volt Operation

Low Voltage from Transformer
for Safety



One Year Warranty on all Parts

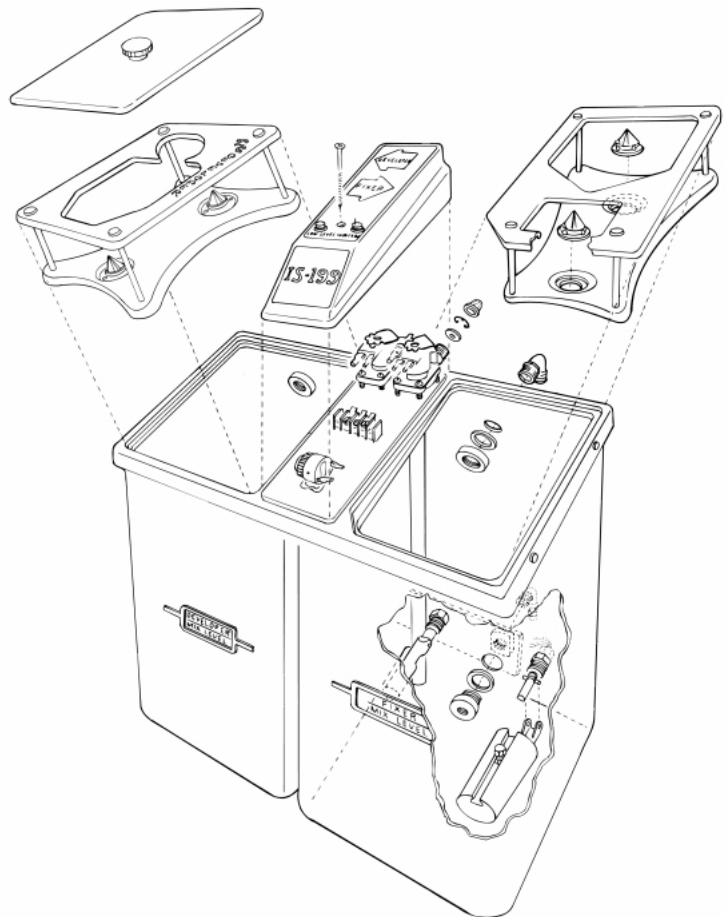

The Versatile Chemical Mixer Mixes the Following Chemicals:

■ Medical X-Ray

■ Industrial X-Ray

■ Graphic Arts

■ Microfilm


The **IS-199** incorporates a unique combination of features which make it the "Mixer of Choice" for both the user and the serviceperson.

The Specific Gravity Mixing System allows the user to start a mix at anytime. The mixer responds to the addition of concentrated chemicals and adds the correct amount of water automatically.

All Electronic Components have been chosen for dependability. Straight-forward circuit designs eliminate the need for circuit boards thereby increasing reliability.

The Low Level Warning System indicates that solution level is 2.0 gallons or less. The warning system resets automatically when the user begins a 5 or 10 gallon mix.

INSTALLATION SITE SPECIFICATIONS

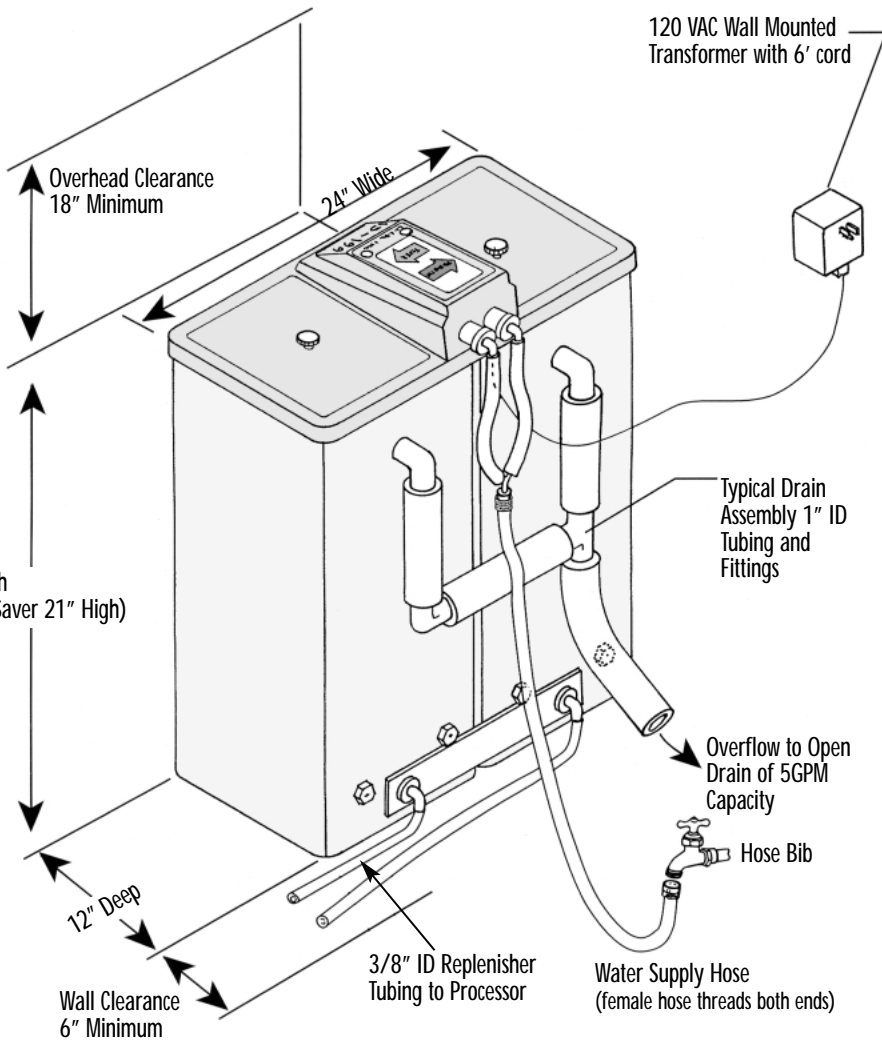

Installer to supply:

- 3/8" Replenisher Tubing
- 1" ID Drain Tubing & Fittings
- High Pressure Water Supply Hose (FGHT Fittings)

TECHNICAL SPECIFICATIONS

Electrical:

120 VAC, 1/2A, 50/60Hz
(240 VAC 1/4 A, 50/60Hz is optional)

Water

18-125 psi, 3GPM or more

Dimensions, Standard Model

30" H x 24"W x 12"D

Dimensions, Space Save

21" H x 24"W x 12"D

Capacity, Standard Model

12.5 Gallons, Dev. & Fixer

Capacity, Space Save

7.5 Gallons, Dev. & Fixer


OTHER INFORMATION

All mixers made to order, allow up to 2 weeks for delivery. Mixers also available in 240 volt internal transformers. Rebuilt mixers available upon request. For more information contact White Mountain Imaging.



ACCESSORIES

Templates for other bottle configurations are available by request on a custom order basis.

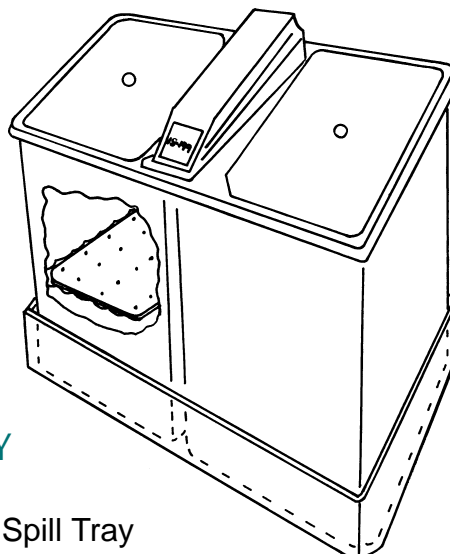


SPILL TRAY

The IS-199 Spill Tray conveniently fits under the IS-199 Chemical Mixer to provide Secondary Containment for potential leaks or spills

FLOATING LID

The IS-199 Floating Lid is designed to extend developer life by reducing aerial oxidation. The IS-199 Floating Lid is especially useful in low volume situations where premature oxidation may occur

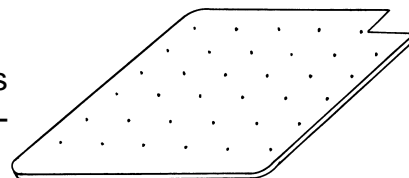

MIXER TEMPLATES

T2 Developer/T2 Fixer
WMI Cine Developer/Fixer
Kodak Developers/Fixers
Repackaged Developers

AGFA Developers/Fixers
Cronex Developer/Fixer
Autex Developer/Fixer
Simon Developer/Fixer

Raytech Diagnostics
<https://raytechdiagnostics.com/>
613-799-4171
sales@raytechdiagnostics.com